

MULTI-REAGENTS HPTLC DERIVATIZER

CHROMAJET
DS20
Reference: BS130.700

- Ensures perfect Repeatability and Accuracy for your TLC analyzes
- Spraying on Precise Zones pre-defined by the operator
- Perfect regularity of the shape of the jet
- Almost total absence of Aerosol formation
- Spray up to 4 reagents on one plate
- Significantly reduces reagent consumption
- Suitable for all common reagents
- Automatic rinsing system
- Stand-alone or PC-controlled, intuitive operation
- User Protection and Security
- Meeting the requirements of GMP/GLP
- For any TLC/HPTLC plates up to 200 x 200 mm



The ChromaJet DS20 is a microprocessor-controlled sprayer for transferring the derivatization reagents onto TLC/HPTLC plates.

The advantages over manual spraying are a significantly reduced reagent quantity, a very uniform spray pattern and the avoidance of aerosol formation.

The application of the reagents is carried out with highest precision and is reproducible at any time.


This device, allowing TLC operations to be carried out in accordance with GLP prescriptions, protects the health of the user against any chemical risk (inhalations and splashes)

The ChromaJet DS20 can operate stand-alone or connected to a PC.

With a special software, different spraying methods can be created, saved and loaded.



Rack for reagents bottles

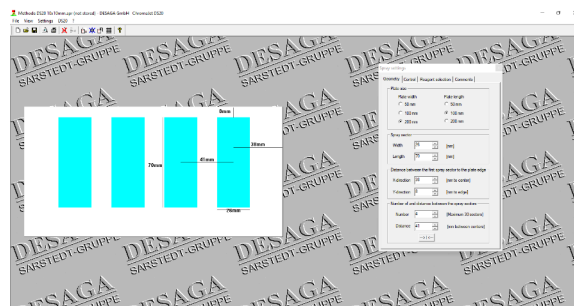


CHROMAJET DS20

MULTI-REAGENTS HPTLC DERIVATIZER

TECHNICAL SPECIFICATIONS

Derivatization area 200 x 200 mm
 Jet speed 50 - 200 mm/s
 Spray cycle 1 – 5
 Spray grid 2 - 4 mm in 0.1 mm steps
 Air throughput 70 - 100 %
 Reagent throughput 50 - 240 s/mL
 Filling rate, rinsing rate 6 - 30 s/mL
 Rinsing operation 1 – 5
 Dimensions (W x H x D) 440 x 405 x 275 mm
 Weight 15 kg

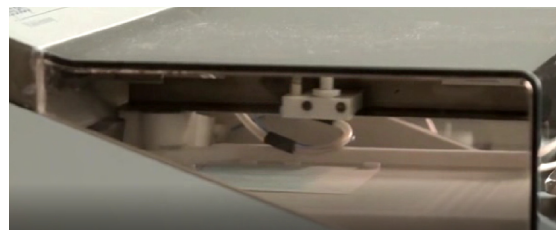

Spraying Nozzle

REFERENCES FOR ORDER

Reference	Description
BS130.700	ChromaJet DS20, reagent spraying device, 230 V, incl. 1 pack of 50 tubes 40 mL, software
BS130.701	ChromaJet DS20, reagent spraying device, 110 V, incl. 1 pack of 50 tubes 40 mL, software
BS130.720	Filter, 1 pack of 20 pcs.
BS130.725	IQ/OQ documents for ChromaJet DS20
BS.63796	Tubes with flat base, screw cap, 40 mL, 1 pack of 50 pcs.
BS130.732	Exhaust hose, viton with connector, 2.50 m

Technical changes reserved. Bionis is not responsible for any misprints, errors which may result in any losses, claims or costs

BIONIS

6 rue de la Prevote
 78550 HOUDAN - France
 Tel: +33 1 75 25 62 05
 E-mail: contact@bionis.fr

DISTRIBUTOR: