

MULTI-WAVELENGTH HPTLC VISUALIZER

HP-UVIS NXG
Reference: BS147.001

- UV light at 254nm (tube), 310 nm (LED),
345 nm (LED) and 366 nm (LED)**
- Low energy consumption thanks to LED lamps**
- Includes UV safety switch**
- Preparative TLC mode**
- 2 inclined plate tables with non-slip surface,
positioned one above the other, equipped with
UV-protection shield**
- 3-sided closed case**
- For TLC/HPTLC plates up to 200 x 200 mm**



It is now possible to examine individual wavelength ranges.

You can choose between 310 nm, 345 nm and 366 nm, or use all three wavelengths together for your test.


The spectrum of the original fluorescent tube can be reproduced using LED technology (all 3 wavelengths active). This allows you to work according to laboratory requirements and subsequently detect the wavelength more accurately.

There is now also a safe preparative mode where you can open the protective screen to 45° and take marks and samples with protective equipment.

If the protective screen is opened more than 45°, the lights switch off immediately.

TECHNICAL SPECIFICATIONS

Detection area:	200 x 200 mm
Dimensions (W x H x D):	325 x 480 x 290 mm
Weight:	11.5 kg
1x UV tube 254 nm	
2x UV LED 310 nm	
1x UV LED 345 nm	
1x UV LED 366 nm	

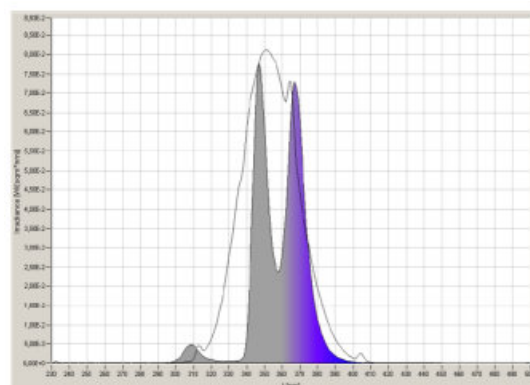

Technical changes reserved. Bionis is not responsible for any misprints, errors which may result in any losses, claims or costs